




compact steam steriliser

instruction manual

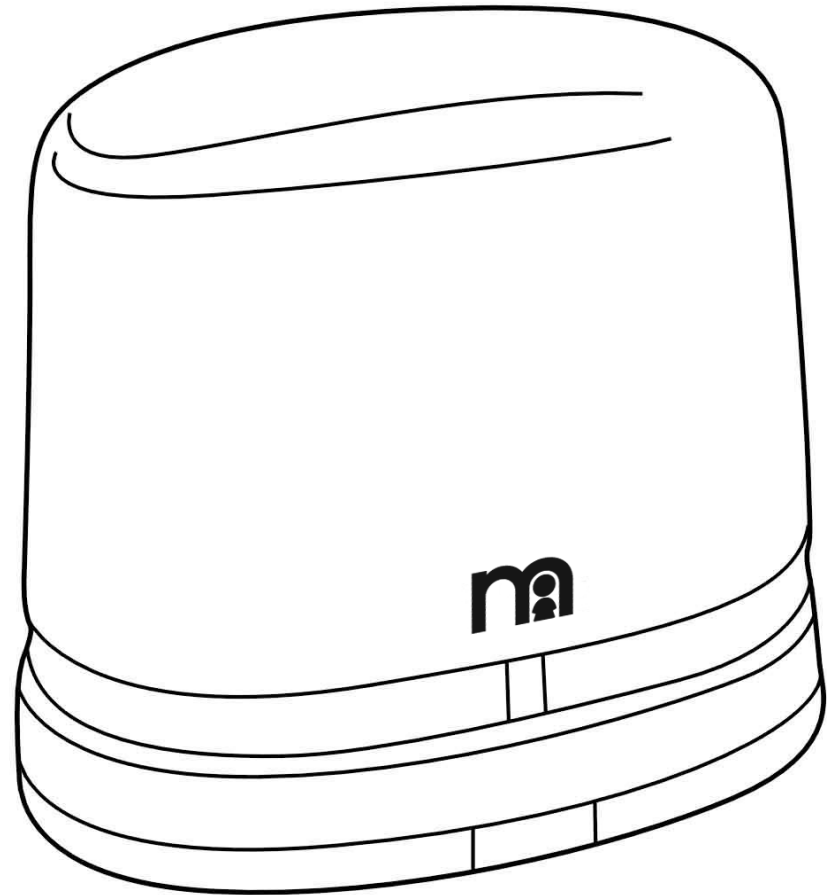

Please keep a note of our company name and address.

Made for Mothercare, Watford, England WD24 6SH

www.mothercare.com

STERILISER INSTRUCTIONS

The Mothercare two bottle Compact Steam Steriliser does not require chemicals for its operation. Sterilisation is achieved by surrounding the bottles and accessories in steam which kills the bacteria.

This Steriliser is designed so that in its compact form, it can be stored conveniently or becomes the ideal unit for travelling. The cable and plug are easily stored in the base for this purpose.

When ready for operation it will sterilise one or two bottles together with teats and accessories in approximately 10 minutes.

LOADING AND OPERATION OF THE STERILISER

Using the measuring cup pour 20ml of water into the tray of the base unit.

NOTE: Do not exceed this volume of water, it is very important that the correct amount of water is used (use exactly 20ml).

Using too little water means there is a risk that the unit will not fully sterilise the contents whereas using too much will prolong the cycle.

Before placing bottles and accessories into the steriliser they should be washed in warm soapy water and rinsed thoroughly with clean warm water.

Place the central section of the steriliser onto the base unit, place all the bottle accessories into the central section and then fit the tray in place. Now place the two bottles upside down onto the tray and replace the upper section, the steriliser is now ready for use.

Plug the unit into a mains socket and the red indicator will illuminate to show the steriliser is connected to the supply.

To start the sterilisation cycle simply press the button on the front of the unit, the green indicator will now illuminate and remain on throughout the cycle, after approximately 10 minutes the green light will go out and the cycle is complete. (please note the unit cannot be switched off once the cycle is completed).

After use wipe the inside surface of the unit clean and dry.

IMPORTANT

The first time you plug your steriliser into the mains supply the green light will illuminate and the steriliser will complete a full cycle. This is quite normal and is due to our testing procedures during manufacture. Please allow the unit to cool prior to reuse.

Do not touch the heater tray during operation. Do not remove cover during operation and when removing it after the sterilisation process is complete allow any condensation to drain back into the body.




CONDENSATION WILL BE PRESENT IN THE UNIT AND THE BOTTLES – THIS IS QUITE NORMAL.

Do not obstruct the steam vents

Always allow the contents to cool before removal.

Always pour off excess water in such a manner that water cannot enter the switch or indicator apertures.

ALWAYS ENSURE THERE IS WATER IN THE TRAY BEFORE COMMENCING STERILISATION PROCESS.

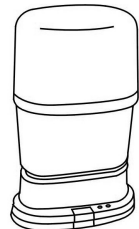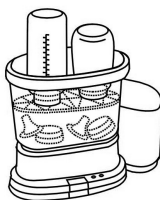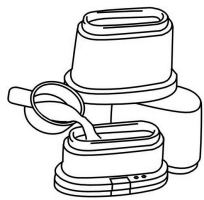

CARE OF YOUR STERILISER

The base unit MUST NOT BE IMMERSSED IN WATER. After unplugging the unit from the mains supply, the Steriliser parts can be wiped clean and dried with a non abrasive cloth.

Scale residue found on the water tray can be removed using a dilute solution of one part vinegar to two parts water. This should be poured into the heater tray and unit turned on without bottles. Repeat two normal cycles with fresh water to rinse, wipe the tray with a damp cloth.

SERVICE

Your Mothercare Compact Steam Steriliser is guaranteed for 24 months from the date of purchase. Should your steriliser fail to operate, check that the fuse in the plug has not blown and that the mains supply is switched on, if problems persist you can call the Mothercare Helpline on 08453 30 40 30.

MAINS PLUG INFORMATION

Your Mothercare Product is supplied with a fitted plug, if for any reason it is to be replaced please follow the procedures below, alternatively ring our helpline on the above number.

If a 13 Amp Plug (BS1363) is used, a 3 Amp fuse must be fitted: or if any other type of plug is used a 5 Amp fuse must be fitted either in the plug or Adaptor or at the distribution board.

When connecting the wires to the plug make sure the blue wire is connected to the terminal marked 'N' or coloured black and connect the brown wire to the terminal marked 'L' or coloured red.

PLEASE NOTE:- Neither wire must be connected to the earth terminal of a 3 Pin Plug.

If the plug supplied is of the non-rewirable (moulded on) type and socket outlets in the home are not suitable for the plug supplied it should be cut off and an appropriate plug fitted as described above.

PLEASE NOTE:- The plug severed from the mains lead must be destroyed as a plug with bared flexible cord is hazardous if engaged in a live socket outlet.

Use a 3 Amp Fuse approved by ASTA to BS1362.

Always replace the Fuse Cover.

Never use the Plug with the Fuse Cover omitted.

Spare Fuse Covers are available from Mothercare Ltd.

The Supply Cord of this appliance cannot be replaced by the user. For all repairs return to our Service Department at the above address.

